

4PT™

4 POLE TECHNOLOGY



Regal 344A


AUSTRALIA

www.meritsaustralia.com

Regal 344A




Regal 344A XL

Features

- Loop handlebars offer comfortable finger or thumb throttle operation
- 4 wheel adjustable suspension system allowing for more rider control and comfort
- Large 14" wheels which provide greater ground clearance
- Larger 4 pole motor and larger controller delivering much more power for climbing steep driveways or hills
- Built-in half speed safety switch which operates automatically when turning to provide extra safety
- New improved, fully adjustable leather-feel seat, allowing maximum rider comfort
- Separate emergency braking system for added security
- With new easy to access 12 volt socket
- Available in wide seat version (Regal 344A XL)

Dimensions cm/in

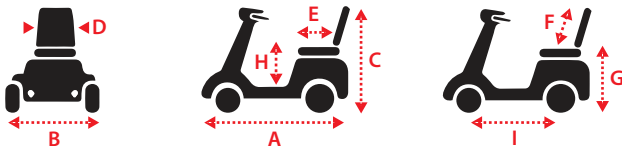

Overall Dimension : A x B x C	152 x 71 x 135cm/60 x 28 x 53"
Seat Dimension : D x E x F	51 x 46 x 51cm/20 x 18 x 20"
Seat Dimension (Regal 344A XL): D x E x F	80 x 46 x 51cm/31.5 x 18 x 20"
Seat Height from Ground : G	56-64cm/22-25"
Seat Height from Deck : H	40-48cm/16-19"
Wheel Base : I	104cm/41"

Weight kg/lb

base	battery/ea	seat	total with battery	capacity
81/ 178	25/55	19/42	150/330	227/500



Technical Specs

Max Speed	12kph/7.5mph
Range	40-50km/26-32mi
Turn Radius	167cm/65.7"
Ground clearance	10cm/4"
Motor	DC 24V, 650W (Rating)/2200W (Max)
Controller	Rhino 160A
Battery	12V/60-80A/H x 2pcs
Charger	8A off-board
Gradient	10°
Front Wheel	14" pneumatic
Rear Wheel	14" pneumatic
Brake	intelligent, regenerative, electromagnetic brakes

merits
AUSTRALIA

www.meritsaustralia.com

Merits Australia is committed to constantly improve and develop products, we reserve the right to alter specifications without prior notice. Pictures for illustration purposes only.